



PRISM MEDICAL Sling Selection

PREMIUM SLINGS

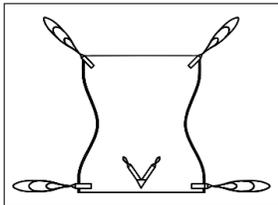


COMFORT RECLINE SLING

The Comfort Recline Sling is designed for the amputee client but it is also suitable for clients who have to remain in the sling for extended periods. The spacer material is fully breathable and stretches to mold to the shape of a chair. The sling features elasticated sections which allow the straps to be folded safely out of the way when the client is seated.

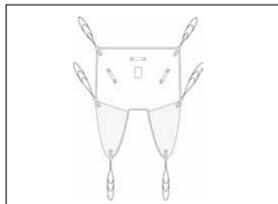


UNIVERSAL SLING

The Universal Sling is the most commonly supplied sling for general transfer purposes. The universal is easy to fit yet is very supportive. It can be used in a divided leg or cross over position.

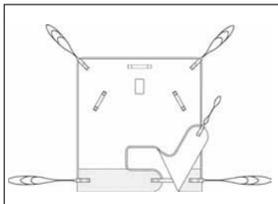

HAMMOCK SLING

The Hammock Sling provides maximum support to the sacral area and is used in cases where the universal sling is unsuitable because the individual has poor body tone. The sling can be used in a divided leg or cross over position; this also provides good support for amputees.

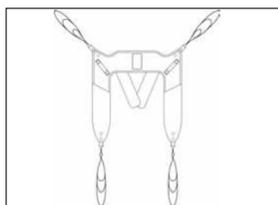

DUAL ACCESS SLING

The Dual Access Sling is designed to be used with either a ceiling lift or a stand aid. It gives access to be able to remove and refit lower garments for toileting transfers while the client is in the sling.



CUSTOM SLING

Although we offer over 100 standard slings, everyone is unique and their needs are different. That is why we will customize any of our standard slings to help suit your needs better.

You can modify an existing sling design, like adding an extra strap, or we can start right from scratch.

Before creating a custom sling, our engineering department reviews the design to ensure that there will be no loss of integrity with the changes.


Once the custom sling has been made by our in-house sling manufacturing, we weight test each custom sling to ensure its strength.

Common customizations include adding: loops, straps and padding.

CLIENT SPECIFIC SLING

Prism Medical is also pleased to offer client specific slings. With Infection Control becoming an area of greater focus and attention, Prism Medical has the solution. Client Specific Slings are designed to be used by one individual and disposed of when; a) the patient/resident no longer requires or b) the sling is soiled.

The Client Specific slings are available in a variety of styles and sizes. These slings are made of a strong, non-woven material that conforms to the body to improve comfort and safety. And because these slings are intended for only one client they help prevent infections from spreading, as well as eliminating laundry costs.

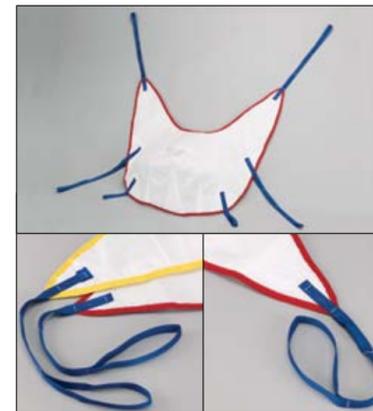

WE'RE THE COMPANY THAT BRINGS YOU A Wide Selection of Slings

EXCEPTIONAL FIT AND DURABILITY



SLINGS TO ACCOMMODATE VIRTUALLY ANY LIFT, POSITIONING AND TRANSFER REQUIREMENT

Prism Medical slings are recognized by health care professionals throughout North America for exceptional fit and durability. Our range of sling styles and sizes will accommodate virtually any lift, positioning and transfer requirement. And, if our standard selection does not provide a sling solution, we will custom design a sling system to your specifications.

Prism Medical slings are manufactured from durable nylon mesh, polyester fabrics and quilted padding to provide contoured support, comfort and safety. As well, these long-wearing, quick-drying slings contain no metal parts, making them convenient to launder and easy to return to service, without the need to locate mislaid inserts.

| Size | Height | Weight |
|----------|---------|-----------------|
| Junior | <4' | <110 lbs |
| Small | 4'-5'6" | 95 lbs-150 lbs |
| Medium | 5'-6' | 125 lbs-250 lbs |
| Large | 6'-7' | 250 lbs-400 lbs |
| X-Large | 6'-7' | 400 lbs-480 lbs |
| XX-Large | >6.5' | >480 lbs |



