

Medcare Low Pro Lifts




Medcare Low Pro Lifts are designed to be used by caregivers to transfer non-weight bearing patients without the strain of manual transfers. The Low Pro's inverted caster design allows it to clear the frame of any manufactured bed or cart while providing smooth travel. This lift is available in two weight capacities including 475 lbs. and 600 lbs. for bariatric patients. For added convenience, an optional scale feature is available on the Low Pro Lift-N-Weigh models to easily weigh a patient simultaneously while transferring.

The expandable base of Medcare's floor lifts permits for lifting around various obstacles such as wheelchairs, beds, tubs, and toilets, creating a safe and user-friendly transfer system.

Item Number:

- 450001: Low Pro Lift 475 - lifts up to 475 lbs.
- 450003: Low Pro Lift-N-Weigh 475 - lifts up to 475 lbs.
- 450019: Low Pro Lift 600 - lifts up to 600 lbs.
- 450013: Low Pro Lift-N-Weigh 600 - lifts up to 600 lbs.

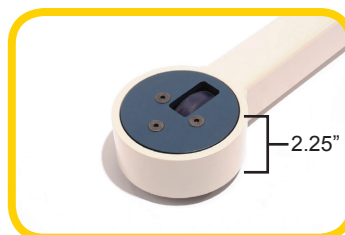

Features:

- *Patented caster design allows for fit under the lowest beds while still minimizing the push/pull force required*
- *Easily transfers patients to and from chairs, wheelchairs, and toilets*
- *Floor and bed transfers can be accomplished with ease*
- *Designed to work with a variety of Medcare slings to "customize" the transfer*

“Lifts are an important part of Medcare's ZBI (Zero Back Injury) System, which reduces costs associated with workers-compensation claims and medical expenses.” ”



Low Pro Lifts (cont.)



Hand Control



Scale Option

Medcare's patented scale option is available on the Care Lift-N-Weigh. This allows for standing transfers to be performed while recording a patient's weight. Whether the patient is standing straight, at an angle, or even holding onto the handles, the scale accurately weighs patients.



Hanger Bar Stopper Option

The Medcare Hanger Bar Stopper prevents sling loops from sliding off of the hanger bar hooks.

Item Number: 400611

Battery Charger



Battery

Technical Specifications

Overall Length:	47"	Height of Wheel Base:	2.25"
Outside Width (legs closed):	27"	Power Source	Rechargeable 24 Volt, 5 Amp Battery
Inside Width (legs open):	36"		

Product Highlights:

- Low profile design allows for access under beds and carts with only 2.25" of ground clearance
- Dual battery design allows stand to be in constant use
- 360° rotation for easy patient positioning
- Arched boom for added head clearance
- Hand control allows caregiver to be positioned next to patient during transfer
- 4-point, non-rigid hanger bar design allows sling to conform to patient for added comfort



Medcare Products

151 East Cliff Road
Burnsville, MN 55337

Toll-free: 800 695 4479
Fax: 952 300 2100

www.MedcareLifts.com

© 2014 Medcare Products Inc.

Rev 031714