



MOTION Bathtub

Unique Relaxing Bathing Experience

Elegant Design Supporting Mobility of Your Residents

MOTION Bathtubs are appropriate for all stages of mobility and compatible with any appropriate stand aid or mobile floor lift. Ergonomic reclining design allows for a safe and comfortable transfer for your resident, the unique tub contour enables optimal sliding from any wheelchair while the rim and shaped seat provide zero obstructions for entry. The MOTION Bathtub features electronic height adjustment providing a low entry height for independent transfers from a wheelchair if needed.




MOTION Bathtub Options



Air-Jet Whirlpool System

Creates a pleasing effect on the skin, promoting a more relaxed tranquil bathing experience.

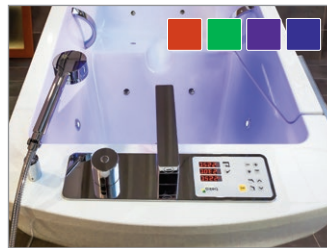

Colour Light Therapy System

Transforms the bathing water into a vibrant colour spectrum creating a stimulating atmosphere.



Disinfection System

Offers increased protection against cross infection with approved in-house cleaning fluids.




Music Sound System

Underwater speakers generate gentle vibrations of music into the water aiding in relaxation.



Automatic Tub Filling

Saves time for the caregiver by auto-filling the tub with warm water for immediate immersion.

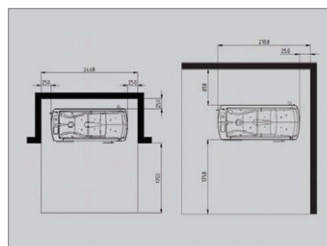

Little Space Required

Saves up to 500 mm (19.7 in) in space compared to standard height adjustable tubs.

Caregiver and Resident Benefits



Sufficient space for the upper body, especially for larger clients. Ergonomical seat and back shapes with hygienic opening for better cleaning access.

Door latch with safety locking for maximum safety.



MOTION is operated with a control panel keypad, with three temperature displays and activation for disinfection.



Extra large door opening facilitates the entry with mobile lifts as well as the self-entry.



Generous foot zone with a pre-filling volume of 98 litres (25.8 gallons) and additionally chromed brass holding grips for safe hold during bathing.



The low entry height of only 490 mm (19.3 in) enables an easy self-entry even from the wheelchair.



Clear view of the room and caregiver because of a no-barrier console.



Motion E is operated with a Cable Remote Control. (Optional with the MOTION)



A transfer with an active or passive lift is possible with the free standing carrier frame.

MOTION Bathtub

Specifications

Safe Working Load (SWL)

Max patient weight: 220 kg (485 lb)

Patient and water: 350 kg (772 lb)

Dimensions (L x W x H)

Tub length, entry position: 1760 mm (69.3 in)

Tub length, tilting position: 1907.8 mm (75.1 in)

Tub width, outside: 860 mm (33.9 in)

Lowest entry height: 490 mm (19.3 in)

Working height: 810 mm / 1010 mm (31.9 in / 39.8 in)

Capacity

Without patient: 270 kg (71.3 gal)

Lifting Range

200 mm (7.9 in)

Inclination Angle

23°

Water Drainage

75 L / minute (19.8 gal / minute)

Filling

Volume, without client: appx. 270 L (71.3 gal)

Pre-filling volume, foot zone: 98 L (25.9 gal)

Filling time: 40 L / minute (15.75 gal / minute)

Power Supply

110V, 60 Hz

Product Weight Without Packaging

175 kg (386 lb)

For additional technical data, please refer to the user manual.

Part Numbers

810150540: MOTION E

810150550: MOTION

Accessories / Options

Description	MOTION E	MOTION
Cable remote control	inclusive	810150555
Automatic tub filling / water stop	810150561	810150560
Music sound system	810150572	810150572
Air-jet whirlpool system	810150571	810150570
Emergency function at power failure	810150587	810150587
Colour light therapy system	810150576	810150575
Disinfection system	810150580	810150585
Second holder for hand shower	910150586	910150586
Comfort neck pillow	810150565	810150565



HANDICARE

485 Millway Avenue, Unit 2, Concord, ON, L4K 3V4

Toll Free: 877.304.5438 | Fax: 416.260.0256

customerservice.canada@handicare.com

www.handicare.ca

10888 Metro Court, St. Louis, MO 63043

Toll Free: 866.891.6502 | Fax: 314.692.7858

customerservice.northamerica@handicare.com

www.handicareusa.com

Are you ready to provide
your residents with a unique
relaxing bathing experience?

To learn more, call 1-888-637-8155 or
visit www.handicare.ca / handicareusa.com